

# Predator® Zirconia



RESTORATIVE

Dentistry and photography courtesy of Dr. Gildo Santos



Cuts through zirconia in as little as 10 seconds.


**Predator® Zirconia** diamond burs are manufactured using a robust bonding technology process that provides maximum cutting efficiency and effortless access through tough zirconia for both endodontic treatment and crown removal. Each bur has a unique black protective coating which insulates the nickle plating ensuring a significant reduction in friction and heat build-up while ensuring superior cutting efficiency and allowing the diamond particles to remain bonded to the bur shank. Often only one **Predator Zirconia** bur is required to access or remove a zirconia-based restoration.



**PREDATOR ZIRCONIA KIT**

2EA tapered fissure, football round, ball and inverted cone burs (medium & fine grit)

## PREDATOR® ZIRCONIA DIAMONDS

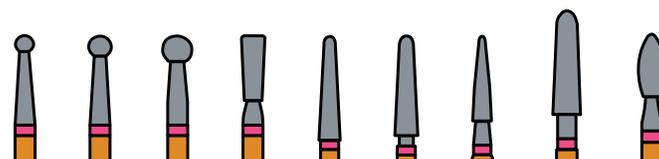

US NO.	801	801	801	807	856	856	862	878K	368
DIAMETER 1/10MM	014	018	023	018	018	022	014	021	023
H/LENGTH (MM)	1.4	1.8	2.3	5.0	9.0	9.0	7.0	8.0	5.0
MEDIUM	●	●	●	●	●	●	●	●	●
FINE	●	●	●	●	●	●	●	●	●

**CLINICIAN'S CHOICE®**  
DENTAL PRODUCTS INC.

[www.clinicianschoice.com](http://www.clinicianschoice.com)