

CounterFit™ II

Multi-Purpose Replication Silicone

Overcoming many clinical problems associated with alginate, such as tearing, voids, lack of impression detail, and single pour restriction, **CounterFit™ II** is a silicone alginate replacement material designed for all alginate-type impressions.

- Flows beautifully to capture exquisite detail for all alginate-type impressions
- Extremely high tear strength compared with traditional alginate and it cures to an ideal durometer to ensure interproximal tags do not fold during reinsertion
- Once polymerized, can be poured more than once and reused, if necessary
- Thixotropic to help prevent patient gagging



CounterFit II: What to Expect

CounterFit II uses the highest quality fillers, resulting in a thixotropic consistency that flows beautifully. Combining the right mixture of catalyst, surfactant and silicone oil, **CounterFit II** produces better surface wetting for exceptional impressions every time.

Sufficient Working Time for Zero Drags, yet a Fast Intraoral Cure Time

Impression drags are caused by having insufficient working time to properly fill a tray and place it intraorally. With **CounterFit II**'s advanced, heat-activated chemistry, working time is not compromised yet once the impression tray is in the mouth, **CounterFit II** reacts to the higher temperature, resulting in a fast setting time, without drags.

Higher Tear Strength and Dimensional Stability

CounterFit II's advanced catalyst system and the proper combination of chain and cross-linking silicones produces higher tear strength, faster elastic recovery and a higher degree of dimensional stability. **CounterFit II** impressions are stable so there is no need for immediate pouring and the impressions can be re-poured more than once, as well as over an extended period of weeks or months. **CounterFit II** impressions will not shrink or swell, ensuring greater dimensional accuracy and detail over alginate impressions for better fitting appliances and restorations. Once set, **CounterFit II**'s extremely smooth surface makes it an ideal temporary matrix material for provisionals that require virtually no finishing or polishing.


Applications

CounterFit II is ideally suited for any application requiring occlusal accuracy.

- Opposing arch impressions
- Bleaching tray & athletic mouthguard fabrication
- Study models
- Orthodontic models, appliances, and retainers
- Matrix material for long span provisional crown and bridge restorations
- TMJ splints
- All alginate-type impressions

MAXIMUM WORKING TIME:	1:10 min.
RECOMMENDED MINIMUM INTRAORAL SET TIME:	1:30 min
TOTAL SET TIME:	2:30 min.

Great flow properties lead to capturing hard and soft tissue detail. Highly accurate articulation using polyvinyl siloxane bite registration enhances the overall quality and accuracy of your appliance restoration.

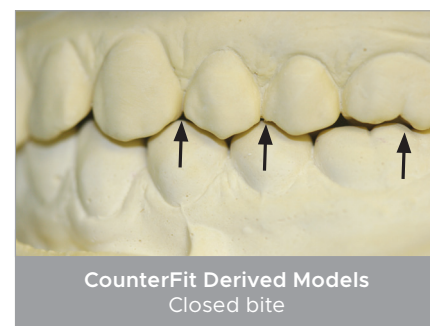

COUNTERFIT II MULTI-PURPOSE REPLICATION SILICONE 50ML 2-PACK (275602)
Contains: 2 x 50mL Cartridges, 6 x 7.5 Mixing Tips, Instructions/SDS

COUNTERFIT II MULTI-PURPOSE REPLICATION SILICONE 50ML 12-PACK (275612)
Contains: 12 x 50mL Cartridges, Instructions/SDS *mixing tips not included

COUNTERFIT II MULTI-PURPOSE REPLICATION SILICONE 50ML 36-PACK (275636)
Contains: 36 x 50mL Cartridges, Instructions/SDS *mixing tips not included

COUNTERFIT II MULTI-PURPOSE REPLICATION SILICONE 380ML MEGAMIX CARTRIDGE (275382)
Contains: 1 x 380mL Cartridge, Instructions/SDS *mixing tips not included

COUNTERFIT II MULTI-PURPOSE REPLICATION SILICONE 120ML 2-PACK (275603)
Contains: 2 x 120mL Cartridges, 6 Dynamic Mixing Tips, Instructions/SDS

NEW