

Detail™

Pre-Impression Cleansing Gel



Detail™ is a 10% EDTA pre-impression cleansing gel that guarantees a clean, decontaminated prep site and, consequently, an accurate impression everytime.

- Lifts organic debris and contaminants like sulfur and methacrylate residue in just 7-10 seconds; leaving behind a squeaky-clean impression surface.
- Lowers the surface tension of the preparation so your impression material contacts the surface more intimately.
- For the best and most efficient impression procedure, use Detail in conjunction with Affinity™ Hydroactive Impression Material.

Fig. 1 – Contaminated preparation site - preventing an accurate and properly cured impression.

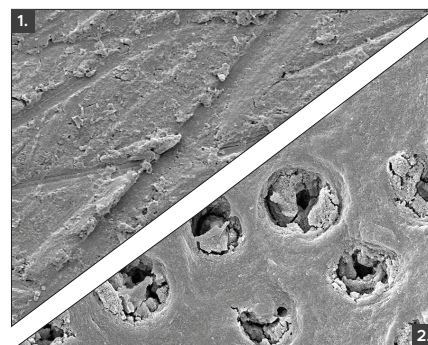



Fig. 2 – Clean preparation site. In just 7-10 seconds, Detail removes contaminants and surface debris, leaving the dentin tubule plugs in place, and allowing for a clean, decontaminated preparation site.

CLINICAL TECHNIQUE

Dentistry courtesy of Dr. Brad Carson. Photography courtesy of Gildo Santos, DDS



Dentistry and photography courtesy of Steve Tsotsos, DDS.

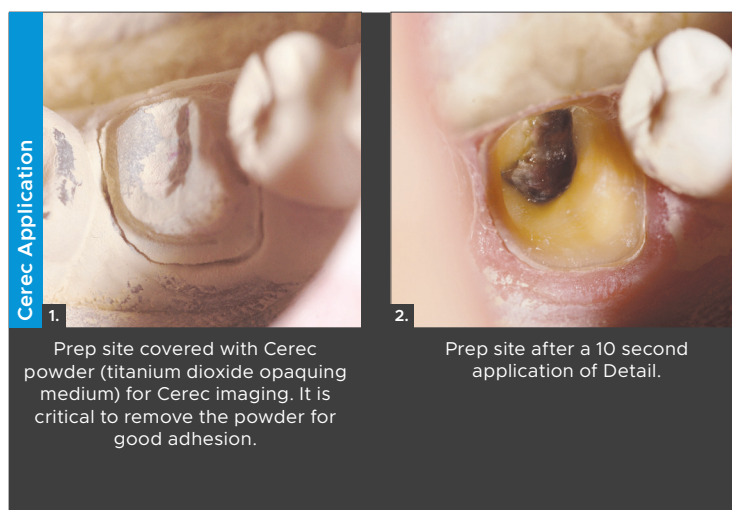

DETAIL PRE-IMPRESSION CLEANSING GEL 4-PACK (035904)

Contains: 4 x 3mL syringes Detail, 50 x flocked needle tips, instructions/SDS



TISSUE GOO/DETAIL TRIAL PACK (642403)

Contains: 2 x 3mL syringes Tissue Goo, 2 x 3mL syringes Detail, 25 x 20 gauge needle tips, instructions/SDS

To view the latest news, education and more, visit:

www.clinicianschoice.com   

 **Clinician's Choice®**
Designed for better dentistry.