

Evanesce™

Nano-Enhanced Universal Restorative



Strong. Simple. Beautiful.

Simple or complex. Single shade or layering technique, anterior or posterior. Whatever your technique or indication, **Evanesce™ Nano-Enhanced Universal Restorative** is designed for truly natural, life-like esthetics.

Exceptional Handling

Evanesce, a nano-hybrid composite material, has a stiffer viscosity, yet its non-sticky formula and ease of manipulation make it highly adaptable. Once placed, Evanesce can be sculpted without slumping to save on finishing time.

A Complete Restorative System

High strength, Evanesce is 82% filled and boasts a high compressive strength, low wear, and low shrinkage to give you confidence for both anterior and posterior restorations.


Evanesce is available in 30 VITA shades and 3 opacities and polishes to an ultra-high, long-term shine in seconds.

Evanesce UX

When a lower viscosity formula is required, try Evanesce UX, now available in 7 popular universal shades. With an 85% universal/body opacity, Evanesce UX shades can be used for single shade restorations or combined with other Evanesce opacities and Enamel FX shades for multi-layer cases. Evanesce UX is available in shades: A1, A2, A3, A3.5, B1, C2, and BL1.


Evanesce FX

To further enhance esthetic enamel characterization, Evanesce is available in three non-VITA "Enamel FX" shades: Enamel White, Enamel Incisal, and Enamel Clear.

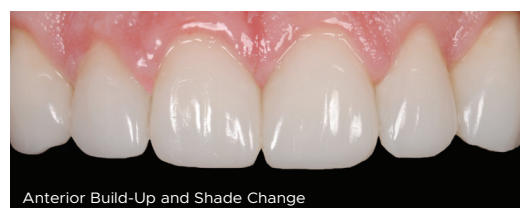

Dentistry and photography courtesy of
Marcos Vargas, BDB, DDS, MS

Evanesce Opacities

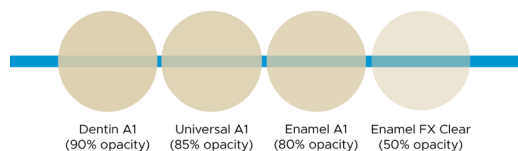

Evanesce FX Enamel White

- 70% opacity
- Reinforced microfill
- Used to create a lingual shelf and to repair small incisal defects



Evanesce FX Enamel Incisal

- 60% opacity
- Nano-hybrid
- Strength and opacity makes it suitable for incisal placement



Evanesce FX Enamel Clear

- 50% opacity
- Reinforced microfill
- Thin clear coat enhances facial lobes



“Developing detailed morphology on anterior restorations is easier with Evanesce due to its no slump handling.”


- Marcos Vargas, BDB, DDS, MS



Request a digital demonstration
for you and your team
findtherestoration.com/evanesce-request-a-demo

Posterior Class I

Dentistry and photography courtesy of Danièle Larose, DDS, AAACD

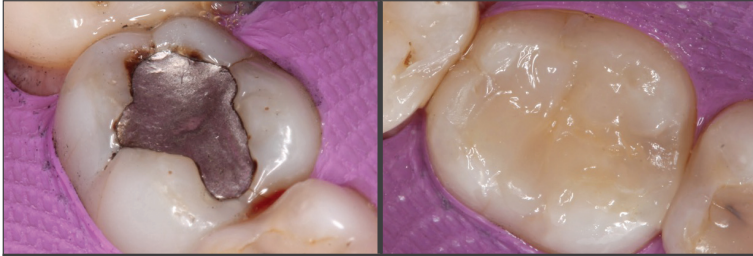

Evanesce packs beautifully into Class I and II preparations when compared to other composites—especially bulk-fill materials. Superior non-slump handling allows you to create occlusal grooves and ridges; minimizing finishing time.

Evanesce Match to Porcelain

Dentistry and photography courtesy of Danièle Larose, DDS, AAACD



Matching composite material to porcelain is a common clinical challenge. The Evanesce range of opacities and outstanding polishability make porcelain matching cases significantly easier.

Multi-Shade Emergency Anterior Case

Dentistry and photography courtesy of Bob Margeas, DDS



Evanesce's putty-like handling allows for easy adaptation during placement and manipulation, yet is absolutely slump-free when you take your instrument or brush away, holding its form until you polymerize.

CLINICAL
Instant
Replay


In a rush?
Watch this quick
technique tip!

POWERED BY
Hands-On
LIVE

Scan to watch

Building Incisal Translucency in Anterior Restorations
David Chan, DDS, AAACD

Evanesce Shades

UNIVERSAL SHADES (85% Opacity)												
A1	A2	A3	A3.5	A4*	B0.5*	B1	B2	C1*	C2	C3*	BL1*	BL3*
UX SHADES** (85% Opacity)												
A1UX*		A2UX*		A3UX*		A3.5UX*			B1UX*		BL1UX*	
DENTIN SHADES (90% Opacity)												
A1*		A2		A3*		A3.5*			B1*		C4*	
ENAMEL SHADES (80% Opacity)												
A1		A2		A3*		A3.5*		B0.5*		B1		B2*
ENAMEL  SHADES												
ENAMEL Incisal				ENAMEL White					ENAMEL Clear			

*Available in single dose only

**UX indicates a slightly lower viscosity



EVANESCE SYRINGE TRIAL KIT (325901)
Contains: 5 x 4g syringes in Universal shades:
A1, A2, A3, A3.5, B1, instructions/SDS

EVANESCE SINGLE DOSE 10pk
Contains: 10 x 0.25g single dose,
instructions/SDS

A4 Universal (322505)
B2 Universal (322507)
C1 Universal (322508)
C3 Universal (322510)
BL1 Universal (322511)
BL3 Universal (322513)
B0.5 Universal (322536)
A1 Enamel (322514)
A2 Enamel (322515)
A3 Enamel (322516)
A3.5 Enamel (322517)
B0.5 Enamel (322537)
B1 Enamel (322518)
B2 Enamel (322519)
A1 Dentin (322520)
A2 Dentin (322521)
A3 Dentin (322522)
A3.5 Dentin (322523)
B1 Dentin (322524)
C4 Dentin (322526)
Enamel FX White (322534)
Enamel FX Incisal (322535)
Enamel FX Clear (322533)

EVANESCE SINGLE DOSE 20pk
Contains: 20 x 0.25g single dose,
instructions/SDS

A1 Universal (322501)
A2 Universal (322502)
A3 Universal (322503)
A3.5 Universal (322504)
B1 Universal (322506)
C2 Universal (322509)
A1UX (323563)
A2UX (323564)
A3UX (323562)
A3.5UX (323565)
B1UX (323567)
BL1UX (323569)

EVANESCE SYRINGES
Contains: 1 x syringe,
instructions/SDS

A1 Universal (321501)
A2 Universal (321502)
A3 Universal (321503)
A3.5 Universal (321504)
B1 Universal (321505)
B2 Universal (321522)
C2 Universal (321506)
A1 Enamel (321507)
A2 Enamel (321508)
B1 Enamel (321511)
A2 Dentin (321517)
Enamel FX Incisal (321515)
Enamel FX White (321514)
Enamel FX Clear (321513)

To view the latest news, education and more, visit:

www.clinicianschoice.com

Clinician's Choice®
Designed for better dentistry.